

Mercury Medical[®]

Your Need ... Our Innovation[®]



Veteran Owned
Small Business

ONLY ONE

FLOWSAFE II[®] +

Disposable BiLevel CPAP System

BiLevel & CPAP Disposable System in One

- **ONLY ONE** O₂ Source
- Less Parts – Easy Set-up

A major leap in product innovation, it takes Patient Care to a whole new level.



www.mercurymed.com

Since 1963

ONLY ONE

FLOWSAFE II[®] +

Disposable BiLevel CPAP System

When lives hang in the balance, you need a disposable BiLevel & CPAP system that lives up to the challenge.

The Flow-Safe II[®]+ Disposable BiLevel & CPAP System from Mercury Medical[®] is the only system that has so many advantages, it gives you an advantage in life-saving emergencies.

In short, the Flow-Safe II+ represents nothing short of a revolution in patient care.



EPAP Pressure Adjustment Dial

INDICATIONS FOR USE:

The Flow-Safe II+ Disposable BiLevel & CPAP device is intended to provide CPAP or BiLevel CPAP pressures to spontaneously breathing patients in the hospital, surgery center and pre-hospital environment.

Clinicians can *now* deliver BiLevel & CPAP therapy with just **ONE** disposable device.



ADVANTAGE 1: BILEVEL/CPAP SWITCH


The Flow-Safe II+ offers the capability of providing both BiLevel or CPAP therapies.

- BiLevel/CPAP switch allows clinicians to provide either therapy mode
- Flow determines CPAP or BiLevel IPAP pressure
- Completely Disposable Single-Patient Use device
- Easy set-up, less parts
- Small compact package easily fits in First-In EMS bags



ADVANTAGE 2: EPAP DIAL

When in BiLevel mode, the EPAP dial allows clinicians to set EPAP pressure

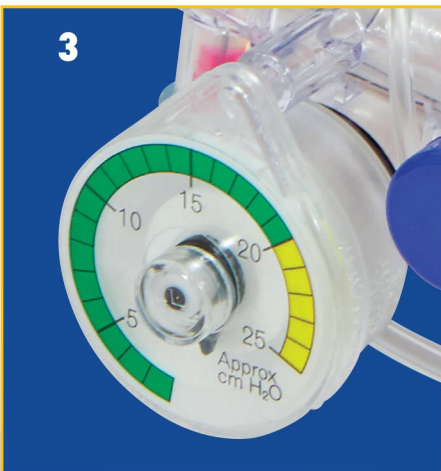

ADVANTAGE 3: BUILT-IN MANOMETER & PRESSURE RELIEF VALVE

The Flow-Safe II+ system features a built-in manometer for verified pressure readings. No assembly of separate apparatus. And the pressure relief safety valve limits pressure to 25 cm H₂O to avoid any unintentional excess pressure.

- Built-in Manometer verifies CPAP pressure
- When providing BiLevel therapy, manometer verifies IPAP and EPAP pressures



ADVANTAGE 4: CPAP SYSTEM CONSUMES 50% LESS OXYGEN THAN ORIGINAL FLOW-SAFE®

The Flow-Safe II+ conserves oxygen while maintaining high FiO₂ delivery. Saves resources and represents a major advantage on long transports. 15 Liters per minute delivers over 10 cm H₂O CPAP and IPAP pressures.



ADVANTAGE 5: ADVANCED MASK DESIGN

The Flow-Safe II+ deluxe face mask with double seal forms a tight anatomical seal without excess pressure, is light and easy to handle. The forehead pads adjust to the patient's anatomy while the stretchy, elastic head harness is easy-to-place with velcro straps that easily adjust for patient comfort.

PRSRT STD
U.S. POSTAGE
PAID
Permit #579
St. Petersburg, FL

Clinicians can *now* deliver BiLevel & CPAP therapy
with just **ONE** disposable device.



FLOWSAFE II[®] +
Disposable BiLevel CPAP System

ORDERING INFORMATION

| Part # | Description | Packaged |
|------------|---|----------|
| #10-574-00 | Flow-Safe II+ BiLevel CPAP with Large Adult Full Face Mask with Straight Swivel Port and Headstrap | 5/Box |
| #10-574-01 | Flow-Safe II+ BiLevel CPAP with Medium/Small Adult Full Face Mask with Straight Swivel Port and Headstrap | 5/Box |
| #10-574-02 | Flow-Safe II+ BiLevel CPAP with Small/Child Full Face Mask with Straight Swivel Port and Headstrap | 5/Box |
| #10-574-03 | Flow-Safe II+ BiLevel CPAP without Mask | 5/Box |

A Revolution in Patient Care

MercuryMedical[®]

11300 - 49th Street North • Clearwater, Florida 33762-4807 USA

ORDERING INFORMATION: 800.835.6633 or 727.573.0088

Fax: 727-573-6040

www.mercurymed.com



Veteran Owned
Small Business

Patents: US 8,522,618 B1 / US 8,893,716 B2 /
US 9,370,635 B2 / Other patent(s) pending