

CLEARVUE[®]

VIDEO LARYNGOSCOPE



CLEARVUE™

The Infinium ClearVue™ VL3R

Reusable Blade Video Laryngoscopes

The Infinium ClearVue™
Video Laryngoscopes are:

COMPACT AND ERGONOMIC

- Lightweight - Unit Is Less Than 350g.
- Quick installation and Removal of Blade
- Quick Shot Camera Button for Recording

LOW COST PER PATIENT


- Reusable Blades Can Be Used Up to 1000 Times

DURABLE

- 2 Years Warranty
- Screen is Reinforced With Anti-Shatter Protection
- Light Weight Long Lasting Rechargeable Lion Battery

HIGH-DEFINITION

- 3.5 Inch Display
- Anti-Fog Lens
- 2.0 Mega Pixel Camera With Video
- Optional HDMI External Monitor

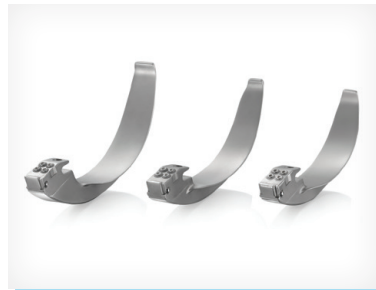

The ClearVue™ VL3R comes standard with size 2, 3, and 4 reusable Mac blades.



Options blade sizes available are: Miller 0, Miller 1, Mac, 1 and Mac 5.

CLEARVUE™

The Infinium ClearVue™ VL3D

Disposable Blade Video Laryngoscopes



The **Infinium ClearVue®** video laryngoscope adopts state-of-the-art image technology and possesses the advantage of a 2.0 megapixel fullview camera with high-resolution monitor. The **ClearVue®** features a rechargeable li-ion battery for extended use. A unique anti-fog capability is also provided. The overall design is compact and ergonomic. A high resolution view of the glottis will enable clinicians first-attempt success while minimizing any chance of a complication during the intubation process.

- Rigid steel chassis construction
- Blade sizes 1, 2, 3, 4, and 5 available
- Vacuum Packed blades
- Low Cost Disposable Blades
- 1 Year Warranty



Size 1, 2, 3, 4, and 5 disposable Mac blades are available for the ClearVue™ VL3D.

CLEARVUE® TECHNICAL SPECIFICATIONS:

DISPLAY

Size	3.5" Full View
Monitor	LCD; Scratch-resistant and anti-glare screen
Resolution Ratio	640*480 RGB
Aspect Ratio	4:3
Video Refresh Rate	30FPS
Illumination	LED

LARYNGOSCOPE BLADE CAMERA INCLUDED

Resolution Ratio	2.0M Pixels
Field Angle	66°
Illuminance	≥800LUX

BATTERY

Type	Rechargeable Lithium Battery
Long working time	≥120minutes
Voltage	3.7V
Capacity	3200mAh
Charging times	More than 300 times
Charging time	4 hours

POWER ADAPTER

Charging port	Micro USB
Input	100-250V,50/60Hz
Output	5V/ 2A

TRANSPORT/ STORAGE

Temperature -10 ~+45

CONDITION

Humidity ≤93
AP 500hPa to 1060hPa

WORKING CONDITION

Temperature 5 ~+40
Humidity 30 ~85
AP 700hPa to 1060hpa

UNIQUE ANTI-FOG FUNCTION

Auto anti-fog function upon powering on

WEIGHT

321 g(disposable system); 376g (reusable system) - light-weight, portable, comfortable

PHOTO/VIDEO FUNCTION

One click image capture 4G internal memory card for video and images USB connection cable for uploading images and video

HANDLE

ABS material - Comfortable ergonomic short handle design