

Ambu® King LTS-D™ Kit

Laryngeal Tube

DATASHEET

The Ambu King LTS-D laryngeal tube is a disposable, simple-to-use alternative airway device that provides superior patient ventilation. The King LTS-D allows the passage of the gastric tube through a separate channel.

The King LTS-D is designed with a straightened, beveled distal tip that assists in directing the airway posterior to the larynx and into the upper esophagus. Due to this unique configuration, there is minimal risk of the device entering the trachea.

A single port inflates both proximal and distal cuffs. The proximal cuff stabilizes the King LTS-D and seals the oropharynx. The distal cuff blocks entry of the esophagus, reducing the possibility of gastric insufflation. The King LTS-D achieves ventilatory seal pressures of 30cm H₂O or more. Multiple distal ventilatory openings and bilateral ventilation eyelets facilitate air flow.



Kit includes one syringe and one pack of sterile lubricating jelly.

KEY BENEFITS

Seals in the esophagus and oropharynx to provide positive pressure ventilation

Single inflation port

Separate gastric access channel

Phthalate and latex free

Color-coded full range of sizes

Ambu

SIZING INFORMATION

Tube Size	Pediatric				Adult		
	Size 0	Size 1	Size 2	Size 2.5	Size 3	Size 4	Size 5
Connector Color	Transparent	White	Green	Orange	Yellow	Red	Purple
Patient Criteria	<5 kg	5-12 kg	12-25 kg 90-115 cm	25-35 kg 105-130 cm	4-5 feet (122-155 cm)	5-6 feet (155-180 cm)	greater than 6 feet (>180 cm)
Recommended Cuff Volume	10 ml	20 ml	35 ml	40-45 ml	50-60 ml	70-80 ml	80-90 ml
Maximum Cuff Pressure	60 cm H ₂ O						
External Diameter of the Tube	9 mm	9 mm	14 mm	14 mm	17.6 mm	17.6 mm	17.6 mm
Bronchoscopy Via Ventilation Lumen	< 3.0 mm	< 3.0 mm	< 4.0 mm	< 4.0 mm	< 6.0 mm	< 6.0 mm	< 6.0 mm
Gastric Tube	10 Fr	10 Fr	16 Fr	16 Fr	18 Fr	18 Fr	18 Fr




King LTS-D Kit Pediatric

ORDERING SPECIFICATIONS

Item No.	Size	Type
KLTS430	0	Ambu® King LTS-D™ Kit with Suction Port, Size 0
KLTS431	1	Ambu® King LTS-D™ Kit with Suction Port, Size 1
KLTS432	2	Ambu® King LTS-D™ Kit with Suction Port, Size 2
KLTS4325	2 ½	Ambu® King LTS-D™ Kit with Suction Port, Size 2.5
KLTS433	3	Ambu® King LTS-D™ Kit with Suction Port, Size 3
KLTS434	4	Ambu® King LTS-D™ Kit with Suction Port, Size 4
KLTS435	5	Ambu® King LTS-D™ Kit with Suction Port, Size 5

Questions? Contact your Ambu representative at **800-262-8462** or visit ambuUSA.com/ambu-king-lts-d

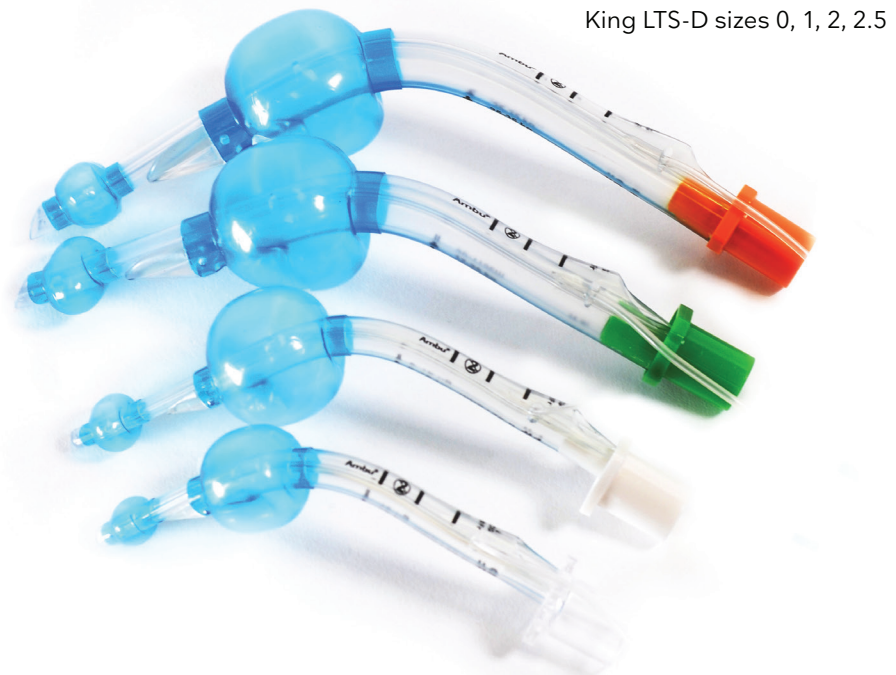

King LTS-D sizes 0, 1, 2, 2.5

Ambu FOREVER FORWARD

Ambu Inc.
6721 Columbia Gateway Drive
Suite 200
Columbia, MD 21046
Tel. 800 262 8462
Fax 800 262 8673
ambuUSA.com

CE 2797 US Rx only