



SAM
MEDICAL

SAM[®] Soft Shell

SHAPE. STRENGTHEN. SUPPORT.

FLEXIBLE COMFORT


SAM[®] Soft Shell Splint is engineered with the same revolutionary C-Curve™ Principle as the original SAM[®] Splint, building strength from curves. With unrivaled flexibility, SAM[®] Soft Shell Splint bends into any simple arc, becoming exponentially stronger and more supportive, immobilizing sprains, strains and minimally displaced fractures.



FEATURED TECHNOLOGY

ADAPTABLE C-CURVE™ PRINCIPLE

SAM[®] Soft Shell Splint leverages the revolutionary C-Curve™ Principle to mold its foam and aluminum layers into structural arcs, strong enough for pre- or post-cast care, while remaining lighter and more transportable than traditional splints that rely on heavy, rigid materials. Bend it back into its original form and it's ready to be remolded.

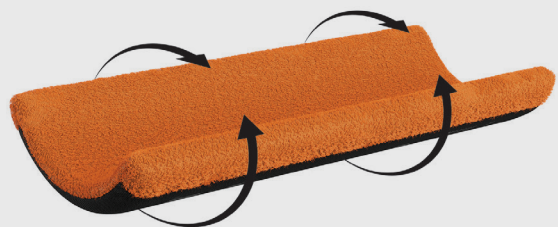




ENGINEERED FOR VERSATILE SUPPORT

- 1 MAXIMUM-COMFORT SURFACE**
Absorbent-brushed terry cloth fabric covers skin-contact surface, providing enhanced comfort.
- 2 MOLDABLE ALUMINUM**
Flat 0-temper aluminum inside the splint's soft outer layer molds easily, allowing for a wide spectrum of applications.
- 3 RADIOLUCENT MATERIAL**
Medical professionals are able to take X-rays and CT-scans without removing the splint from the fractured or injured limb.

- 4 SIZE THAT FITS**
Full product line covers most sizes needed, including XL and sugar tong.
- 5 RESISTANCE CONTROL EDGE**
Firm and extra-firm ends allow different levels of resistance to motion, according to needs. (Not available on the Sugar Tongs.)

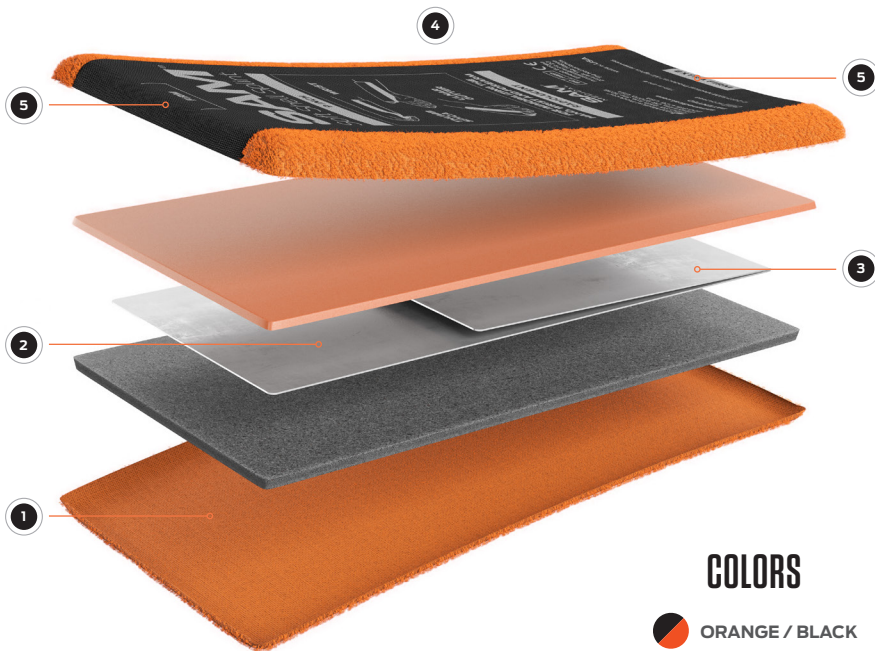

- COLORS**
- ORANGE / BLACK
 - BLACK / BLACK

MODELS



STANDARD
4.25 in wide



XL
5.5 in wide



SUGAR TONG
Pediatric, Adult, Adult XL

PRODUCT SPECIFICATIONS

MODELS AVAILABLE*

9" (SS540-EN)	15" (SS544-EN)
9" XL (SS541-EN)	15" XL (SS545-EN)
12" (SS542-EN)	
12" XL (SS543-EN)	

SUGAR TONG MODELS

Pediatric (ST561-EN)
Adult (ST560-EN)
Adult XL (ST562-EN)

GLOBAL SALES

U.S. Toll Free: 800.818.4726
 U.S. Local: 503.639.5474
 U.S. Fax: 503.639.5425
sammedical.com

Ergo Europe, Prinsessegracht 20, 2514 AP, The Hague, The Netherlands		SAM Medical Products® 12200 SW Tualatin Rd., Ste 200 Tualatin, OR 97062 USA	

* Standard width 4.25 in wide, XL 5.5" wide.
 SAM® is a registered trademark in the United States and other countries.
 ©2019 SAM Medical Products. All rights reserved.