

GRIP-ET

Endotracheal Tube Holder

TUBE HOLDER IMPORTANCE:

Research shows that 5-15% of ET Tubes placed in the pre-hospital setting are either displaced or extubated upon admission to the emergency department. AHA guidelines now recommend the use of commercially available tube holders.

Grip-ET is the solution. Once you've confirmed tube placement with a PosiTube, lock it into place with the only tube holder that leaves 50% of the oral cavity accessible. Unaffected by body fluids and comfortable for patients, the Grip-ET is the ideal choice for stabilizing ET tubes.

MODEL SELECTION:

GT-100XL - 4.5-14mm w/ out Bite Block

GT-202XL - 4.5-14mm w/ Bite Block


KEY FEATURES:

LOCKING-CLIP slides into a patented gear rack which quickly secures tube.

VELCRO STRAP quickly and easily wraps around the patient's head to secure Grip-ET to face.

OFFSET POSITION provides wide access to the mouth. Visualize, suction, or perform oral care without moving the Grip-ET.

RED TWIST LOCK provides the extra security on tube position. The lock pivots to hold the clip in place.

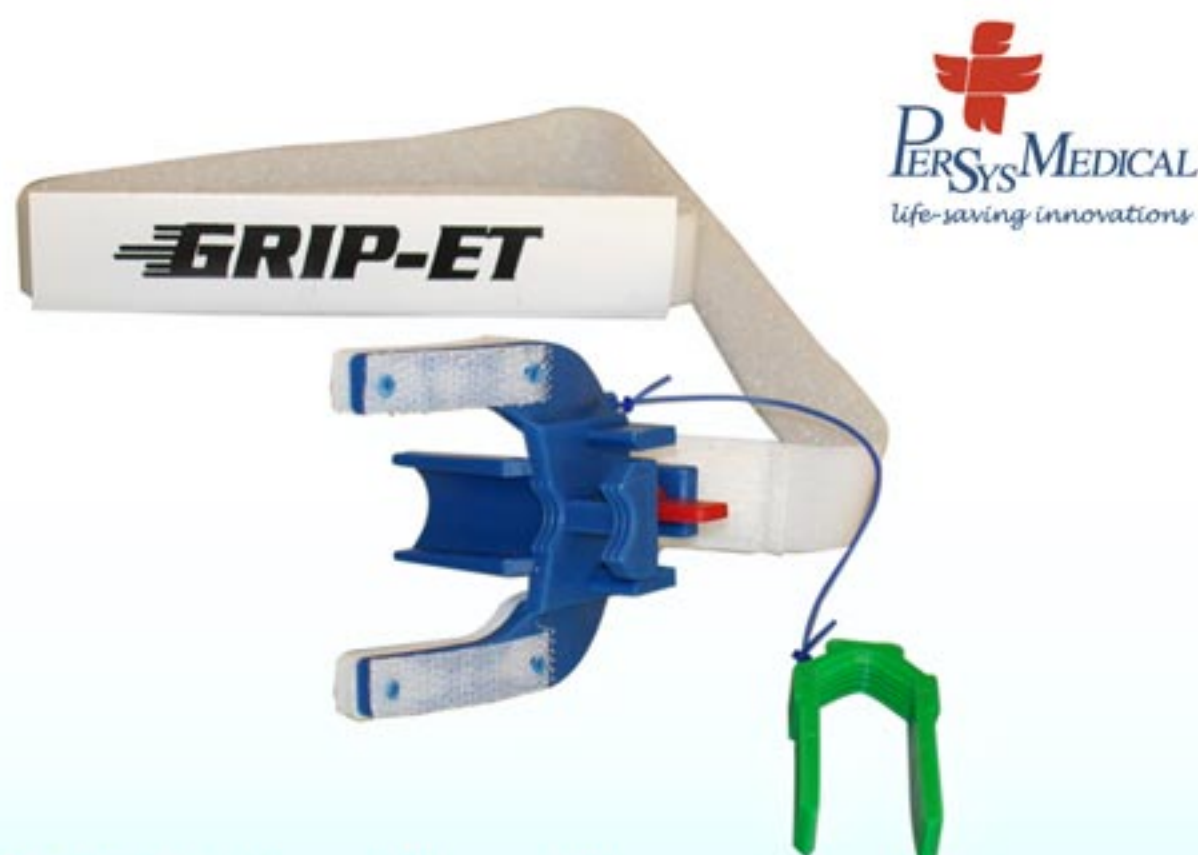

Grip-ET is pictured above with:

- Optional bite block
- Cardboard stick to feed Velcro strap behind neck
- Extra long clips to hold 4.5-14mm tube sizes

BENEFITS:

- Secures any tube 4.5-14mm
- Fits any size LMA, Cobra, King, and Combitube
- Prevents accidental tube movement and extubation
- Available with or without a bite block
- Meets all AHA guidelines
- Secures at either side of the mouth
- Leaves exceptional access to mouth
- May be used with a surgical cricothyrotomy