

# SSCOR DuCanto Catheter




**“A simple and effective tool to support routine and emergency airway management.”**

- James C. DuCanto, M.D.

- Larger internal diameter
- Designed to fit the Oropharyngeal Anatomy
- Can be used in the SALAD technique

The **SSCOR DuCanto Catheter™** is a single-use, non-sterile, rigid suction catheter specifically designed to be used by emergency care providers to aspirate liquids or semi-solids from a patient's upper airway.

It has a **larger inside diameter** than a standard Yankauer suction tip to facilitate the removal of fluids as well as solid material. The larger suction lumen makes it significantly less likely to clog in an emergency situation.

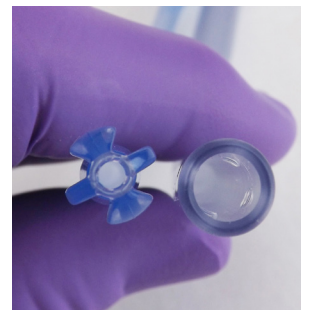

**SSCOR** DuCanto  
Catheter

*Clearing The Airway Is Our #1 Priority*

The **hyper curved catheter** is shaped more like a patient's airway making it easy to position the catheter and work with other complementary devices used to clear a patient's airway.



**SSCOR.com**

The **SSCOR DuCanto Catheter™** is constructed of a rigid PVC material. The larger internal diameter allows for larger pieces of emesis to pass while the lack of thumb port makes suction immediately available when the catheter is connected to the suction source. By removing the vacuum control, the caregiver does not need to locate and keep their thumb correctly positioned over the thumb port to clear a patient's airway.



**SALAD**

Suction Assisted Laryngoscopy  
& Airway Decontamination

Instrumental in the **NEW SALAD Technique.**

The **SALAD (Suction Assisted Laryngoscopy and Airway Decontamination)** technique calls for constant suction provided by a catheter in the upper airway, continuously removing vomit, blood and other materials, around which the practitioner intubates the patient via either direct or video laryngoscopy.

Find more information on **SALAD**  
at [www.sscor.com/salad](http://www.sscor.com/salad).

**SSCOR DuCanto Reference/Part Number: 200-00002**

**Specifications:**

**Shaft diameter: .26" x .55" (6.6mm x 13.9mm)**

**Shaft length: 9.3" (236.22mm)**

For best results, use **SSCOR #43200 OR #43203 patient connecting tubing**. Both have a large internal diameter similar to the **SSCOR DuCanto Catheter**.

Dr. James C. DuCanto is an anesthesiologist with 20 years' experience in private practice who possesses substantial experience in the clinical practice and teaching of difficult airway management, as well as the use of patient care simulators.

SSCOR, Inc.  
11064 Randall St., Sun Valley, CA 91352  
Email: [sales@sscor.com](mailto:sales@sscor.com) or [techsupport@sscor.com](mailto:techsupport@sscor.com)  
Phone: (800) 434 5211 Domestic +1 818 504 4054 International  
[www.sscor.com](http://www.sscor.com) & [shop.sscor.com](http://shop.sscor.com)

**SSCOR** | **DuCanto**  
**Catheter**

*Clearing The Airway Is Our #1 Priority*



**#43200 - 6 Feet**



**#43203 - 3 Feet**

Dr. DuCanto works at the Aurora St. Luke's Medical Center in Milwaukee, Wisconsin where he serves as a staff anesthesiologist, Director of the Simulation Center, and director of the Anesthesiology rotation for students and residents.

**SSCOR**

Made in the USA of U.S. and Imported parts.  
Lit 200-00002 Rev B