

Products designed with a purpose...
...to help you help them.



SURFLO® and SurFlash® I.V. Catheters

An ultra-sharp double-bevel introducer needle plus a medical grade lubricant allow for easier penetration and smoother travel through tissue. Preview chamber in clear sure-grip hub gives immediate indication of proper placement. Flexible catheter with large inner diameter assures good blood flow. All this plus safe, easy-to-handle, tamper-evident packaging make Terumo's I.V. catheters a great clinical value.

LATEX FREE


TERUMO®

TERUMO MEDICAL PRODUCTS

SurFlash: Flashback Feature

With SurFlash, visual confirmation of catheter placement occurs due to a unique flashback feature. Flashback appears between the specially grooved needle and the catheter, providing medical staff with a higher certainty of proper catheter placement. The specially grooved needle also improves the ability to visualize flashback in patients where venipuncture is more difficult. And because SurFlash is made of polyurethane, which is firm at room temperature and softer at body temperatures, patient comfort is enhanced, without compromising proper technique.

NOTE: The needle insertion angle in the illustrations differs from actual insertion angle. When using SurFlash for venipuncture, follow the same procedures and take the same safety precautions as when using Terumo's conventional IV catheter models.

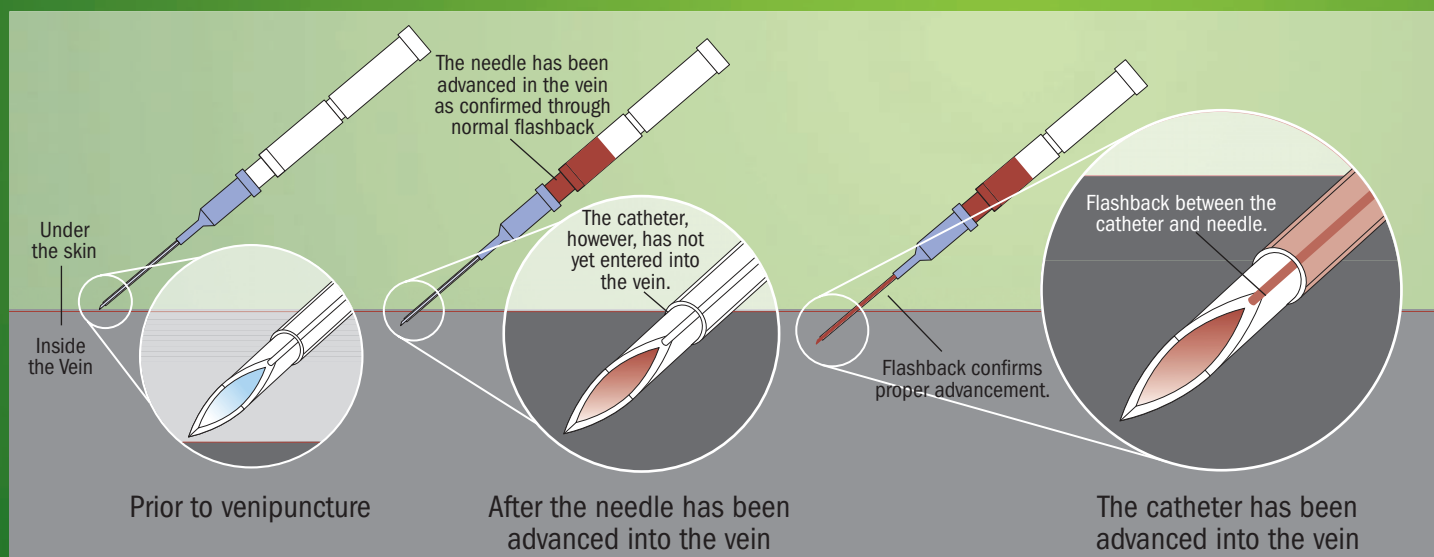

SurFlash Polyurethane I.V. Catheters

Packed: 50 per box, 4 boxes per case (200).

Product Code	Description
SR*FF1451	14G x 2" (Orange)
SR*FF1651	16G x 2" (Gray)
SR*FF1664	16G x 2 1/2" (Gray)
SR*FF1832	18G x 1 1/4" (Green)
SR*FF1851	18G x 2" (Green)
SR*FF1864	18G x 2 1/2" (Green)
SR*FF2025	20G x 1" (Pink)
SR*FF2032	20G x 1 1/4" (Pink)
SR*FF2051	20G x 2" (Pink)
SR*FF2225	22G x 1" (Blue)
SR*FF2419	24G x 3/4" (Yellow)

SURFLO ETFE I.V. Catheters

Packed: 50 per box, 4 boxes per case (200).

Product Code	Description
SR-OX1451CA	14G x 2" (Orange)
SR-OX1632CA	16G x 1 1/4" (Gray)
SR-OX1651CA	16G x 2" (Gray)
SR-OX1664CA	16G x 2 1/2" (Gray)
SR-OX1832CA	18G x 1 1/4" (Green)
SR-OX1851CA	18G x 2" (Green)
SR-OX1864CA	18G x 2 1/2" (Green)
SR-OX2025CA	20G x 1" (Pink)
SR-OX2032CA	20G x 1 1/4" (Pink)
SR-OX2051CA	20G x 2" (Pink)
SR-OX2225CA	22G x 1" (Blue)
SR-OX2419CA	24G x 3/4" (Yellow)

TerumoTMP.com/IVCatheters

US Office

2101 Cottontail Lane
Somerset, NJ 08873
Phone: 732.302.4900
Fax: 732.302.3093
Toll free: 800.283.7866
Fax for Orders: 800.411.5870

Canada Office

Terumo Medical Corporation
C/O UPS - Supply Chain Solutions
1165 Creditstone Road
Vaughan, Ontario L4K 4N7
Phone: 800.668.5489
Fax: 905.760.1215

Ⓢ and TERUMO are trademarks owned by Terumo Corporation, Tokyo, Japan, and they are registered with the U.S. Patent & Trademark Office. SURFLO and SurFlash are trademarks owned by Terumo Medical Corporation, Elkton, MD, U.S.A.

09/07/2012-2404