

# VBM Pelvic Sling

The Pelvic Sling helps to fix and stabilize fractures of the pelvis at the accident scene. Fractures of the pelvis are very painful and the bleeding resulting thereof might lead to a life-threatening condition for the patient.

The integrated inflatable cuff exerts a lateral pressure to the pelvis for best adjusting and stabilization result. This helps to reduce life-threatening inner bleeding and pain significantly.




## SLIDE SPINT / KNEE FIXATION

To slide the sling underneath the patients elvis.  
After removal used as knee fixation (see photo).

## CUFFS

Both lateral cuffs are connected and can be inflated with the hand inflator for an ideal stabilization of the pelvis.

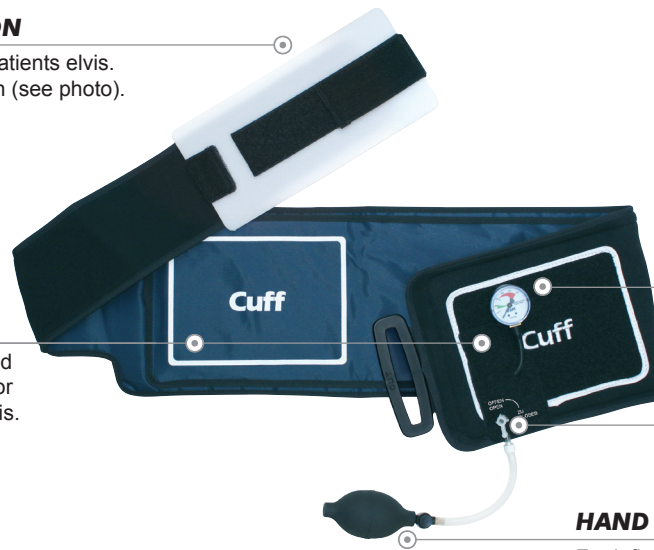

## MRI-PROOF MANOMETER

For indication of the recommended pressure.

## STOP COCK

For a leak free transport.

## HAND INFLATOR

For inflation and deflation of both cuffs.

## ORDER INFORMATION

### Pelvic Sling

with manometer and hand inflator,  
MRI-compatible, for single use

Size	Patient	Pelvis circumference	REF
S	junior	70 – 90 cm	REF 34-10-006
M	standard	90 – 110 cm	REF 34-10-004
L	large	110 – 140 cm	REF 34-10-005

### Pelvic Sling

with manometer and hand inflator,  
reusable

Pelvis circumference	REF
90 – 130 cm	REF 34-10-003

